

DATA SHEET

TG-2022 SUPERBRAVA

SUPERBRAVA SERIES

ELECTRIC EQUIPMENT



Aluminium



Electric



320 Kg
unassisted



Water
resistant



R10




CERTIFIED
20G



- ▷ Control of the stretcher from the operator's position.
- ▷ Load capacity up to 320 kg without assistance.
- ▷ UTAC 20G certification.
- ▷ Manual engine release by means of an engine clutch system.
- ▷ A single battery with extra capacity reaching 50 cycles tested with 100Kg.
- ▷ Secondary control panel allowing automatic unfolding of the legs.
- ▷ Adjustable handrails with 6 positions and capacity to accommodate bariatric patients.
- ▷ Electronic articulated headrest area.
- ▷ Optional bariatric kit with shoulder support for increased patient comfort.
- ▷ Telescopic headrest profile for manoeuvring in reduced spaces.
- ▷ Emergency lights for safety.
- ▷ Integrated carrying handle for front transfers.
- ▷ Electromagnetic compatibility according to R10 standard.
- ▷ Specific model available for incubators.
- ▷ Waterproof and easy-to-clean design.

PRODUCT SPECIFICATIONS

		TG-2022
H	Max. height (mm)	1220
HC	Max. adjustable height (mm)	1090
H2	Min. height (mm)	490
L	Length (mm)	2095
L2	Length closed (mm)	1733
L3	Length open (mm)	2223
W	Width (mm)	590
	Load capacity (kg)	320
	Weight (kg)	78

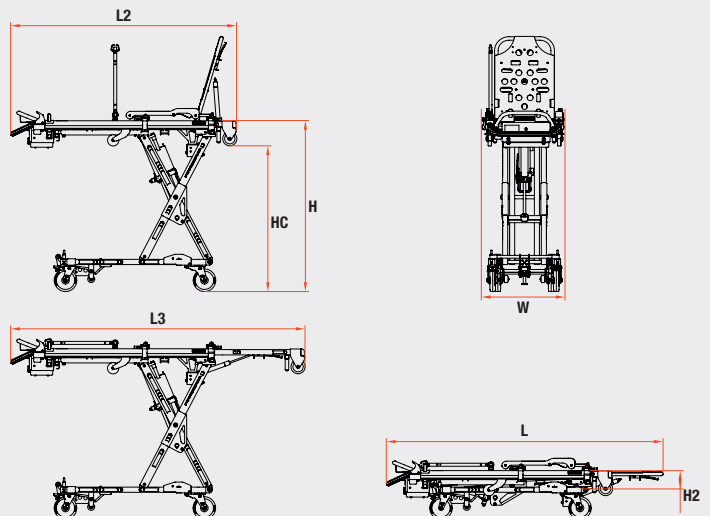

ACCESSORIES



UTILITY TRAY




FIXING FOR CORPULS



PUSHING POLES



BATTERY CHARGER

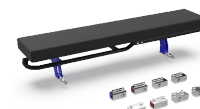

BIARIATRIC KIT



PEDIATRIC ADAPTED
MATTRESS



KARTSANA

Life-moving innovation

Narcís Monturiol, 34
08192 · Sant Quirze del Vallès
Barcelona (Spain)

T.: +34 93 715 86 72
info@kartsana.com
www.kartsana.com